

M is for Mobile EEG Solution

The mobile and portable EEG lab for hospitals and clinics

The Hospital and Clinic Mobile EEG Lab

Ideal for Routine or LTM studies. Continuous or Intermittently Monitored Video EEG made easy with this portable all-in-one system.

We've revolutionized the traditional cart-based hospital EEG system and made it functional for multiple use-cases. Trackit M is designed for rapid deployment in any part of the hospital or clinic setting. Now EEGs can be performed rapidly, anywhere, anytime.

Portable EEG

The entire system fits in a custom case but is flexible enough for use on a wheeled cart, or mounted in a patient room. Technologists can bring EEG to the patient anywhere in the hospital to expand services without large infrastructure costs. Streaming study data can be accessed by physicians anywhere, anytime, via a web browser on any device using our Rendr software.



Rapid Study Setup

The Trackit M includes everything needed in one package for extra mobility:



Infrared HD remote control pan-tilt-zoom video camera



Photic Stimulator



Built-in EEG amplifier



Hand-held event marker



Standard 10/20 Headbox





Rendr Platform cloud software for access anytime, anywhere

rendr[®]
PLATFORM

Start Recording Studies in a Few Clicks with Rendr Acquisition

Capture high-quality data and video while streaming to the Rendr Portal. Our simplified workflow ensures a study can be started with an anonymous patient easily. Intuitive step-by-step instructions on our user-friendly interface provide everything EEG technologists need.

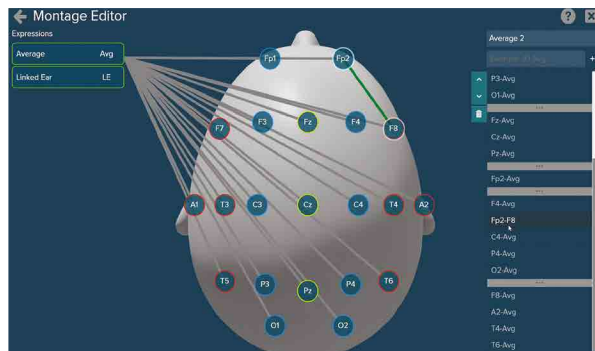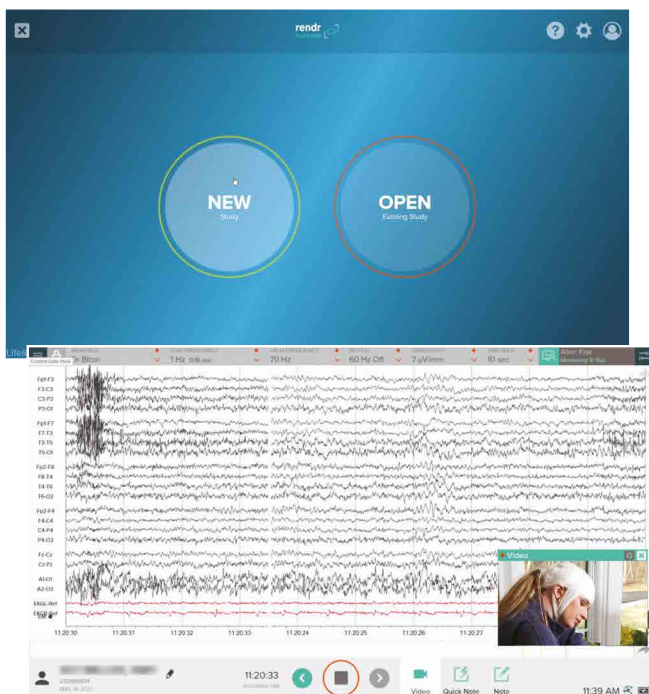


EEG Software Customizability

EEG technologists can create their own color themes, allowing the highest visibility of waveforms and events that are easy on the eyes. Channels can be hidden to highlight only the most critical information.

Create montages depending on what you need. Adjust filters, time-base, sensitivity, and more with our modern control interface.

Notes and events can be added with a keystroke. The duration of these events and notes can be added for focused review at a later time. The pan-tilt-zoom camera can be controlled so the patient is always visible.

Our acquisition software was designed with EEG technologist feedback. Our cloud review tools continue to improve based on what our customers need thanks to our user-feedback collection feature built into the software.

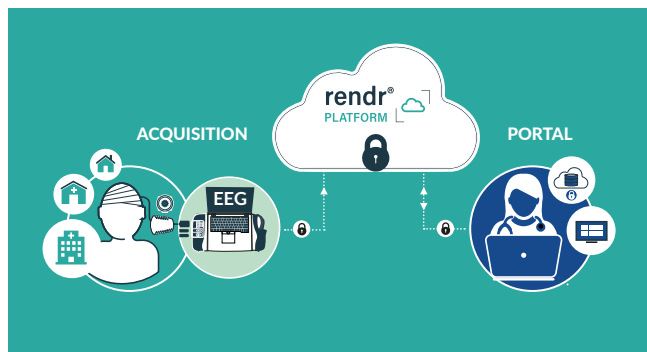


Never Lose a Study

We understand the importance of never losing a study. Recordings are saved simultaneously to a physical storage device and streamed to the cloud so you'll always have a backup. In the event of an issue, our cloud-connected software reports status and alerts for technical problems to our best-in-class support team in real-time. We may know you have an issue before you do, and we have the troubleshooting information at our fingertips to help get you back up and running, fast.

Streaming Data Review, from Anywhere at Anytime

Access and review studies from most internet-connected browsers on laptops, Macs and tablets. Reduce friction and delays for physicians, with portable, on-the-go review. Quickly enter and exit studies while allowing physicians to add notes, documents, and complete their reports - increasing patient throughput.

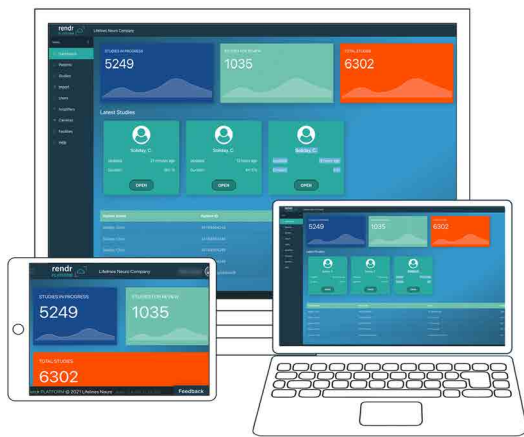

Simplified, Secure Access

The Rendr Portal Practice Management Dashboard allows instant insight into completed or in-progress studies. Increase efficiency by simply sending secure link invitations to verified users who need access to the Rendr Portal. Assign certain roles to users to provide the exact amount of access required for a HIPAA-ready workflow.

Data Protection

Amazon Web Services (AWS), the host of the Rendr Platform, includes HIPAA-ready architecture. AWS includes load sensing and auto-scaling to meet increased bandwidth demand so that as your patient load grows your ability to record and review their studies will too.



Safety

- IEC 60601 compliant for medical electrical equipment
- Two levels of built-in patient protection mean no need to use an isolation transformer for power into the system - can operate off any standard wall AC outlet

World Class Customer Support

Lifelines Neuro is known for the exceptional customer support provided by our EEG-knowledgeable staff, many of whom are registered EEG technologists. Support is available 24/7/365.

EEG Anywhere Subscription

Our flexible subscription plan features the Trackit M bundled with our Rendr Platform software for one fixed monthly fee. With unlimited seat licensing and rapid replacement if your system needs repair, our EEG Anywhere Subscription provides predictable monthly billing, with no large up front capital expenditure, to rapidly scale your business. As your partner, we'll work with you to scale and expand the number of patients you serve.

Turnkey Services

Lifelines Neuro partners with a variety of EEG service providers who can provide turnkey services for hospitals and health systems to expand services.

Contact us today to learn more about expanding services throughout the hospital - ICU, ED, EMU and more with minimal upfront capital investment.